

Assessment findings

in common peripheral vestibular disorders



Touch the future

with Micromedical VisualEyes™
VNG by Interacoustics



Touch
activated
software

User-friendly features

- Remote operation
- Voice and audio cues
- Automatic fixation light
- Seamless video playback
- Customizable protocols



Complete VNG solution for balance assessment

Videonystagmography provides ideal conditions for the observation, measurement and analysis of eye movements. It saves time, reduces costs and avoids the problems of myogenic and electrical noise.

- Modern touch enabled software
- Situational camera for any test and patient interviews
- Customisable user-defined protocols in compliance with BSA standards
- Large eyes on second monitor
- Voice timer with audio cues
- Patient de-identification and PDF reports

interacoustics.com
micromedical.com
academy.interacoustics.com (training materials)


Interacoustics

Assesment findings

Page
















Inferior Vestibular Neuritis (acute)	4
Inferior Vestibular Neuritis (compensated)	5
Superior Vestibular Neuritis (acute)	6
Superior Vestibular Neuritis (compensated)	7
Labyrinthitis (acute)	8
Labyrinthitis (compensated)	9
Meniere's Disease (early stage)	10
Meniere's Disease (late stage)	11
Superior semicircular canal dehiscence	12
Posterior Canal BPPV	13
Horizontal Canal BPPV	14
Anterior Canal BPPV	15
Bilateral Vestibular Loss	16
Vestibular Schwannoma (early stage)	17
Vestibular Schwannoma (late stage)	18

This resource is a tool based on the needs of medical professionals and students that allows quick access to the typical assessment findings in a range of common vestibular disorders. The resource was developed to provide fast, easy-to-use, and always available information which can aid in reaching the correct diagnosis. The information contained within is provided as an information resource only, and should not be used as a substitute for professional diagnosis and management.



Assessment findings
















Inferior Vestibular Neuritis (acute)

Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Normal / not changed	<ul style="list-style-type: none"> Nystagmus present Direction fixed with fast phases towards the good ear Suppresses or vanishes with fixation 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus following Alexander's law Direction fixed Suppresses or vanishes with fixation Linear slow phases 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits Evidence of spontaneous nystagmus may be present on tracings 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits Evidence of spontaneous nystagmus may be present on tracings
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	cVEMP 
Gain: enhanced when the stimulus moves in the same direction of the spontaneous fast phases	Negative	Negative	<ul style="list-style-type: none"> Total response: within normal limits Unilateral Weakness: < 25% Directional Preponderance: <30% Fixation index: <50% 	<ul style="list-style-type: none"> Amplitude: reduced on affected ear or absent cVEMP Latency: within normal limits Ratio: >36%
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Amplitude: within normal limits Latency: within normal limits Ratio: <33% 	<ul style="list-style-type: none"> Gain >0.7 Asymmetry <7% Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain < typically reduced outside of the normative range Phase: Increased phase lead above normal range Symmetry: Asymmetry towards affect ear 	<ul style="list-style-type: none"> Time constant: <10 seconds Gain < typically reduced outside of the normative range Symmetry: Asymmetry towards affect ear

Assessment findings

Inferior Vestibular Neuritis (compensated)


















Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Normal / not changed	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	cVEMP 
Gain: within normal limits	Negative	Negative	<ul style="list-style-type: none"> Total response: within normal limits Unilateral Weakness: < 25% Directional Preponderance: < 30% Fixation index: < 50% 	<ul style="list-style-type: none"> Ratio: >36% or <36% if resolved Amplitude: Reduced amplitude on affected ear or resolved Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: <33% Amplitude: within normal limits Latency: within normal limits 	<ul style="list-style-type: none"> Gain >0.7 Asymmetry <7% Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain within normal limits Phase : within normal limits Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: > 10 seconds but < 60 seconds Gain within normal limits Symmetry: within normal limits

Assessment findings

Superior Vestibular Neuritis (acute)


















Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Normal / not changed	<ul style="list-style-type: none"> Nystagmus present Direction fixed with fast phases towards the good ear Suppresses or vanishes with fixation 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus following Alexander's law Direction fixed Suppresses or vanishes with fixation Linear slow phases 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits Evidence of spontaneous nystagmus may be present on tracings 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits Evidence of spontaneous nystagmus may be present on tracings
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	oVEMP 
Gain: enhanced when the stimulus moves in the same direction of the spontaneous fast phases	Negative	Negative	<ul style="list-style-type: none"> Total response: reduced Unilateral Weakness: > 25% towards affected side Directional Preponderance: >30% towards unaffected side Fixation index: <50% 	<ul style="list-style-type: none"> Ratio: <36% Amplitude: within normal limits Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: >33% Amplitude: reduced on affected ear or absent oVEMP Latency: within normal limits 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain < typically reduced outside of the normative range Phase : Increased phase lead above normal range Symmetry: Asymmetry towards affect ear 	<ul style="list-style-type: none"> Time constant: <10 seconds Gain < typically reduced outside of the normative range Symmetry: Asymmetry towards affect ear



Assessment findings
















Superior Vestibular Neuritis (compensated)

Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Normal / not changed	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	oVEMP 
Gain: within normal limits	Negative	Negative	<ul style="list-style-type: none"> Total response: reduced Unilateral Weakness: > 25% towards affected side or resolved Directional Preponderance: <30% Fixation index: <50% 	<ul style="list-style-type: none"> Ratio: <36% Amplitude: within normal limits Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: >33% or <36% if resolved Amplitude: Reduced amplitude on affected ear or resolved Latency: within normal limits 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain within normal limits Phase : within normal limits Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: > 10 seconds but < 60 seconds Gain within normal limits Symmetry: within normal limits



Assessment findings











Labyrinthitis (acute)

Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
<p>Hearing loss associated with onset of vertigo</p>	<ul style="list-style-type: none"> Nystagmus present Direction fixed with fast phases towards the good ear Suppresses or vanishes with fixation 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus following Alexander's law Direction fixed Suppresses or vanishes with fixation Linear slow phases 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits Evidence of spontaneous nystagmus may be present on tracings 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits Evidence of spontaneous nystagmus may be present on tracings
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	cVEMP 
<p>Gain: enhanced when the stimulus moves in the same direction of the spontaneous fast phases</p>	<p>Negative</p>	<p>Negative</p>	<ul style="list-style-type: none"> Total response: reduced Unilateral Weakness: > 25% towards affected side Directional Preponderance: >30% towards unaffected side Fixation index: <50% 	<ul style="list-style-type: none"> Ratio: >36% Amplitude: reduced on affected ear or absent cVEMP Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: >33% Amplitude: reduced on affected ear or absent oVEMP Latency: within normal limits 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain < typically reduced outside of the normative range Phase: Increased phase lead above normal range Symmetry: Asymmetry towards affect ear 	<ul style="list-style-type: none"> Time constant: <10 seconds Gain < typically reduced outside of the normative range Symmetry: Asymmetry towards affect ear

Assessment findings

Labyrinthitis (compensated)


















Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Hearing loss resolved either partially or completely	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	oVEMP 
Gain: within normal limits Symmetry: within normal limits	Negative	Negative	<ul style="list-style-type: none"> Total response: reduced Unilateral Weakness: > 25% towards affected side or resolved Directional Preponderance: <30% Fixation index: <50% 	<ul style="list-style-type: none"> Ratio: >36% or <36% if resolved Amplitude: Reduced amplitude on affected ear or resolved Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: >33% or <36% if resolved Amplitude: Reduced amplitude on affected ear or resolved Latency: within normal limits 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain within normal limits Phase : within normal limits Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: > 10 seconds but < 60 seconds Gain within normal limits Symmetry: within normal limits

Assessment findings

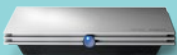Meniere's Disease (early stage)


















Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Characterized by a episodic hearing loss affecting the low frequencies	<ul style="list-style-type: none"> • Not present during remission stage • Present if seen during an attack • Nystagmus can reverse immediately after an attack 	<ul style="list-style-type: none"> • Not present during remission stage • Present if seen during an attack 	<ul style="list-style-type: none"> • Latency: within normal limits • Accuracy: within normal limits • Velocity: within normal limits 	<ul style="list-style-type: none"> • Gain: within normal limits • Symmetry: within normal limits
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	cVEMP 
Gain: within normal limits	Negative	Negative	<ul style="list-style-type: none"> • Total response: reduced • Unilateral Weakness: > 25% towards affected side or <25% • Directional Preponderance: <30% • Fixation index: <50% 	<ul style="list-style-type: none"> • Ratio: >36% • Amplitude: reduced on affected ear or absent cVEMP • Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> • Ratio: >36% • Amplitude: reduced on affected ear or absent oVEMP • Latency: within normal limits 	<ul style="list-style-type: none"> • Gain >0.7 maybe higher than 1.2 • Asymmetry <7% or >7% if increased gain • Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> • Gain >0.7 sometimes > 1.2 • Asymmetry <7% or >7% if increased gain • Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> • Gain within normal limits • Phase : within normal limits • Symmetry: within normal limits 	<ul style="list-style-type: none"> • Time constant: > 10 seconds but < 60 seconds • Gain within normal limits • Symmetry: within normal limits



Assessment findings
















Meniere's Disease (late stage)

Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Persistent hearing loss affecting the low frequencies	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	cVEMP 
Gain: within normal limits	Negative	Negative	<ul style="list-style-type: none"> Total response: reduced Unilateral Weakness: > 25% towards affected side Directional Preponderance: <30% within normal limits Fixation index: <50% 	<ul style="list-style-type: none"> Ratio: ----- Amplitude: cVEMP absent Latency: -----
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: ----- Amplitude: oVEMP absent Latency: ----- 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain <0.7 on affected side Asymmetry >7% towards affected side Catch up saccades: present >50% 	<ul style="list-style-type: none"> Gain < typically reduced outside of the normative range Phase : Increased phase lead above normal range Symmetry: Asymmetry towards affect ear 	<ul style="list-style-type: none"> Time constant: <10 seconds Gain < typically reduced outside of the normative range Symmetry: Asymmetry towards affect ear

Assessment findings

Superior semicircular canal dehiscence


















Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
<p>A 'false' air-bone gap which is greater in the low frequencies</p>	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities May be induced if subject to intense noise or pressure change 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities May be induced if subject to loud noise or pressure change 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	oVEMP 
<p>Gain: within normal limits</p>	<p>Negative</p>	<p>Negative</p>	<ul style="list-style-type: none"> Total response: within normal limits Unilateral Weakness: < 25% Directional Preponderance: < 30% Fixation index: < 50% 	<ul style="list-style-type: none"> Ratio: < 36% Amplitude: within normal limits Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: < 33% Amplitude: within normal limits Latency: within normal limits 	<ul style="list-style-type: none"> Gain > 0.7 Asymmetry < 7% Catch up saccades: not present or < 50% 	<ul style="list-style-type: none"> Gain > 0.7 Asymmetry < 7% Catch up saccades: not present or < 50% 	<ul style="list-style-type: none"> Gain within normal limits Phase : within normal limits Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: > 10 seconds but < 60 seconds Gain within normal limits Symmetry: within normal limits



Assessment findings

Posterior Canal BPPV
















Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Normal / not changed	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities. 	<ul style="list-style-type: none"> No nystagmus present 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	oVEMP 
Gain: within normal limits	Torsional up beating nystagmus towards the under most ear and may reverse in return-to-sitting (RTS). Down beating for apo-geotropic (non-ampullary) PC and rotation away from affected ear.**	Negative	<ul style="list-style-type: none"> Total response: within normal limits Unilateral Weakness: < 25% Directional Preponderance: <30% Fixation index: <50% 	<ul style="list-style-type: none"> Ratio: <36% Amplitude: within normal limits Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: <33% Amplitude: within normal limits Latency: within normal limits 	<ul style="list-style-type: none"> Gain >0.7 Asymmetry <7% Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> Gain >0.7 Asymmetry <7% Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> Gain within normal limits Phase : within normal limits Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: > 10 seconds but < 60 seconds Gain within normal limits Symmetry: within normal limits

**Note: In apogetropic posterior canalithiasis (non-ampullary arm), if the rotary component is present, the nystagmus will be away from the affected ear (e.g., right apo-PC = clockwise and left apo-PC = counter clockwise).
 Modified by 360 Neuro Health (w/permission from Interacoustics).



Assessment findings

Horizontal Canal BPPV

Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Normal / not changed	<ul style="list-style-type: none"> Nystagmus absent or pseudo spontaneous nystagmus may be present 	<ul style="list-style-type: none"> None or pseudo spontaneous nystagmus may be present 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits Evidence of spontaneous nystagmus may be present on tracings 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits Evidence of spontaneous nystagmus may be present on tracings
Optokinetic 	Dix Hallpike 	Head Roll** 	Caloric Test 	cVEMP 
Gain: within normal limits	May evoke horizontal nystagmus	Geotropic or Ageotropic nystagmus present. Geotropic: affected side is side w/stronger intensity of nystagmus. Ageotropic: affected side is weaker intensity nystagmus.**	<ul style="list-style-type: none"> Total response: within normal limits Unilateral Weakness: < 25% Directional Preponderance: <30% Fixation index: <50% 	<ul style="list-style-type: none"> Ratio: <36% Amplitude: within normal limits Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: <33% Amplitude: within normal limits Latency: within normal limits 	<ul style="list-style-type: none"> Gain >0.7 Asymmetry <7% Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> Gain >0.7 Asymmetry <7% Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> Gain within normal limits Phase : within normal limits Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: > 10 seconds but < 60 seconds Gain within normal limits Symmetry: within normal limits
















**Note: In Horizontal Canal (HC) BPPV, the fast phase of the nystagmus beats towards the affected ear -- regardless whether it is geotropic or apogeotropic. Transient geotropic indicates otoconia in the posterior arm of the HC while apogeotropic indicates canalithiasis in the anterior arm or associated with cupulolithiasis. Consider light or heavy cupula or central if nystagmus doesn't fatigue (geo = light cupula; apogeo = heavy cupula).

Modified by 360 Neuro Health with permission by Interacoustics.

Assessment findings

Anterior Canal BPPV












Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Normal / not changed	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits Evidence of spontaneous nystagmus may be present on tracings
Optokinetic 	Dix Hallpike** 	Head Roll 	Caloric Test 	oVEMP 
Gain: within normal limits	Transient down beating nystagmus towards the under most ear. (Rotary component present 50% of the time. Nystagmus doesn't reverse in the RTS)	Negative	<ul style="list-style-type: none"> Total response: within normal limits Unilateral Weakness: < 25% Directional Preponderance: <30% Fixation index: <50% 	<ul style="list-style-type: none"> Ratio: <36% Amplitude: within normal limits Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: <33% Amplitude: within normal limits Latency: within normal limits 	<ul style="list-style-type: none"> Gain >0.7 Asymmetry <7% Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> Gain >0.7 Asymmetry <7% Catch up saccades: not present or <50% 	<ul style="list-style-type: none"> Gain within normal limits Phase : within normal limits Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: > 10 seconds but < 60 seconds Gain within normal limits Symmetry: within normal limits

**Note: The Straight Head Hanging Test (SHHT) is the most reliable test for anterior canal (AC) BPPV. If the rotary component is present, the nystagmus will be towards the affected ear (e.g., right AC = counter clockwise and left AC = clockwise).

Assessment findings

Bilateral Vestibular Loss


















Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Normal / not changed	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits Evidence of spontaneous nystagmus may be present on tracings
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	cVEMP 
Gain: within normal limits	Negative	Negative	<ul style="list-style-type: none"> Total response: < 22 degrees per second Ignore all other parameters 	<ul style="list-style-type: none"> Ratio: ----- Amplitude: cVEMP absent Latency: -----
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: ----- Amplitude: oVEMP absent Latency: ----- 	<ul style="list-style-type: none"> Gain <0.7 bilaterally Asymmetry either >7% or <7% Catch up saccades: present >50% in each direction 	<ul style="list-style-type: none"> Gain <0.7 bilaterally Asymmetry either >7% or <7% Catch up saccades: present >50% in each direction 	<ul style="list-style-type: none"> Gain < reduced outside of the normative range Phase : Increased phase lead above normal range Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: <10 seconds Gain < typically reduced outside of the normative range Symmetry: none

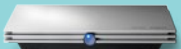

Assessment findings
















Vestibular Schwannoma (early stage)

Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Unilateral sensorineural hearing loss	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits Evidence of spontaneous nystagmus may be present on tracings
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	cVEMP 
Gain: within normal limits	Negative	Negative	<ul style="list-style-type: none"> Total response: within normal limits Unilateral Weakness: < 25% Directional Preponderance: < 30% Fixation index: < 50% 	<ul style="list-style-type: none"> Ratio: > 36% Amplitude: reduced on affected ear or absent cVEMP Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: > 33% or < 36% if resolved Amplitude: Reduced amplitude on affected ear or resolved Latency: within normal limits 	<ul style="list-style-type: none"> Gain < 0.7 on affected side Asymmetry > 7% towards affected side Catch up saccades: present > 50% 	<ul style="list-style-type: none"> Gain > 0.7 Asymmetry < 7% Catch up saccades: not present or < 50% 	<ul style="list-style-type: none"> Gain within normal limits Phase : within normal limits Symmetry: within normal limits 	<ul style="list-style-type: none"> Time constant: > 10 seconds but < 60 seconds Gain within normal limits Symmetry: within normal limits

Assessment findings

Vestibular Schwannoma (late stage)



Pure Tone Audiometry 	Spontaneous Nystagmus 	Gaze Testing 	Saccades 	Smooth Pursuit 
Progressive unilateral sensorineural hearing loss	<ul style="list-style-type: none"> Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Gaze Evoked Nystagmus absent or present only at subclinical slow phase velocities 	<ul style="list-style-type: none"> Latency: within normal limits Accuracy: within normal limits Velocity: within normal limits 	<ul style="list-style-type: none"> Gain: within normal limits Symmetry: within normal limits Evidence of spontaneous nystagmus may be present on tracings
Optokinetic 	Dix Hallpike 	Head Roll 	Caloric Test 	cVEMP 
Gain: within normal limits	Negative	Negative	<ul style="list-style-type: none"> Total response: within normal limits Unilateral Weakness: < 25% Directional Preponderance: < 30% Fixation index: < 50% 	<ul style="list-style-type: none"> Ratio: >36% Amplitude: reduced on affected ear or absent cVEMP Latency: within normal limits
oVEMP 	vHIT lateral 	vHIT vertical 	SHA 	Step Test 
<ul style="list-style-type: none"> Ratio: >36% Amplitude: reduced on affected ear or absent oVEMP Latency: within normal limits 	<ul style="list-style-type: none"> Gain < 0.7 on affected side Asymmetry > 7% towards affected side Catch up saccades: present > 50% 	<ul style="list-style-type: none"> Gain < 0.7 on affected side Asymmetry > 7% towards affected side Catch up saccades: present > 50% 	<ul style="list-style-type: none"> Gain < typically reduced outside of the normative range Phase: Increased phase lead above normal range Symmetry: Asymmetry towards affect ear 	<ul style="list-style-type: none"> Time constant: < 10 seconds Gain < typically reduced outside of the normative range Symmetry: Asymmetry towards affect ear

References

- Agrawal, Y., Carey, J P., Della Santina, C.C., Schubert, M C., Minor, L B. (2009). Disorders of balance and vestibular function in US adults: data from the National Health and Nutrition Examination Survey, 2001-2004. *Archives of Internal Medicine*, 25, 169 (10): 938-944
- Baloh, R W. (2003) Vestibular Neuritis. *The New England Journal of Medicine*, 348: 1027-32
- British Association of Otorhinolaryngologists Head and Neck Surgeons (2002). Clinical effectiveness guidelines: acoustic neuroma. BAO-HNS Document 5
- Minor, L B. (2005) Clinical manifestations of superior semi circular canal dehiscence. *Laryngoscope*, 115: 1717-1727
- Parnes, L S., Agrawal, S K., & Atlas, J (2003). Diagnosis and management of benign paroxysmal positional vertigo (BPPV). *Canadian Medical Association Journal*, 169 (7) :681-93
- Salte, A N. & Plontke, S K. (2010). Endolymphatic hydrops: pathophysiology and experimental models. *Otolaryngologic Clinics of North America*, 43(5): 971-983
- Thompson, T L., & Amedee, R. (2009). Vertigo: a review of common peripheral and central vestibular disorders. *The Ochsner Journal*, 9 (1) 20-26



Interacoustics

Interacoustics is a world leading provider of diagnostic solutions in the field of hearing and balance assessment. Since 1967 we have designed and manufactured our innovative diagnostic solutions for the world of audiology with a constant focus on providing our customers with quality, dependable products.

With you at all times

We operate in over 100 countries worldwide through a carefully selected network of distributors and service centres. Purchasing an Interacoustics product guarantees not only a quality product, but also direct access to our highly professional training and support service.

Learn more at:
interacoustics.com/academy



Interacoustics A/S

Audiometer Allé 1
5500 Middelfart
Denmark

T +45 6371 3555
F +45 6371 3522

info@interacoustics.com
interacoustics.com

